

DX One Series



DX One series is built with modular construction and wide range of features to increase patient's comfort and care giver's efficiency. It has been engineered with the most modern manufacturing processes and the result is the combination of efficient performance along with durability which makes it the easiest choice of the user.



*The product images shown are for illustrative purpose only and not be an exact representation of the product. Accessories shown in the product image may not be part of the standard configuration.

We are interested in growing relationship



Technical Specifications

Overall Dimension*	2145 mm x 1020 mm (L x W)
Platform Frame*	2080 mm x 955 mm (L x W)
Mattress Platform*	1900 mm x 860 mm (L x W)
Mattress Platform Height* (Without Mattress)	535 mm to 780 mm
Backrest**	0° - 70°
Kneerest**	0° - 45°
Safe Working Load	200 Kgs
Patient Load bearing capacity^	175 Kgs


Electrical Parameters

Power In (Mains Power Supply)	100 - 240 V ~, 50/60 Hz., Max 5A
Electrical Shock Protection	Class I
Degree of Shock Protection	Type B
Liquid Ingress Protection	IPX6
Duty Cycle	10 % (Two minutes for every eighteen minutes)

Variants

Features	DX4311 380	DX4311 360	DX4211 380	DX4211 360
Knee rest	Motorized	Motorized	Manual	Manual
Trendelenburg**	0° - 10°	0° - 10°	0° - 9°	0° - 9°
Reverse Trendelenburg**	0° - 7°	0° - 7°	0° - 9°	0° - 9°
Auto Contour	Motorized	Motorized	-	-
Coverage span	Full length	¾ th length [§]	Full length	¾ th length [§]
Angle indicator for backrest, upper leg, trendelenburg/reverse trendelenburg	Available	-	Available	-

Add-On Accessories

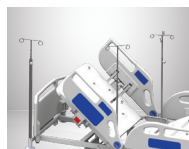

Heavy Duty and
Telescopic IV Pole



O2 Cylinder Cage



Modular
Chart Holder



Orthopedic
Traction Pulley



Heavy Duty
SS Lifting Pole



3rd Set of Safety
Side Railings

- All variants are available with X-Ray permeable backrest (XR) feature as an add-on accessory, which are factory fitted.
- Battery backup[#]

* ±10 mm Engineering variation. ** ± 5° Engineering variation. ^When bed is configured with standard accessories. #Factory fitted.

[§] Can be converted into full length by adding 3rd set railing.

Note: Option of fixed polymer moulded handle available in lieu of removable handle.

Midmark (India) Private Limited

Head Office: Art Guild House, FOF, B-16, "B" Wing, 4th floor, Phoenix Market City, LBS Road, Kurla (West), Mumbai - 400 070
Toll Free: 1800 22 8020 | Tel: +91 22 4915 3000 | Fax: +91 22 4915 3100 | www.midmark.in | india@midmark.com

DX One Series: 09:08.01.24

Standard Features

- Motorized backrest & height adjustment
- Polymer moulded split type safety side rails
- Polymer moulded head & foot panels, easily removable without using any tools
- Dual side manually operated CPR lever
- Four corner safety rubber buffer of 100 mm diameter
- 125 mm diameter single wheel high grade synthetic body castors with diagonal locking mechanism
- Inbuilt SMPS
- 100 mm thick four section mattress
- Removable Handle
- Urine bag holder
- SS heavy duty IV pole (two hooks)

Safety Standards

General safety for medical beds	IEC 60601-2-52
General safety for medical electrical equipment	IEC 60601-1-1
EMI/EMC	IEC 60601-1-2

Accessories and color shown in the picture may not be part of standard configuration. Technical alterations reserved in view of our continuous improvements. Midmark India Pvt. Ltd. is the licensed user of the trademark "Midmark India" and the logo. Reproduction of this catalogue in any form is liable for legal action.