

ICU Solutions

EX6000 Plus



EX6000 Plus
with full length coverage

A bed with state-of-the-art technology and an attractive appearance, EX6000 Plus offers quality beyond comparison, blended with wide range of features such as patient comfort, caregiver convenience, ease of mobility and is engineered to international standards. This new bed from Midmark India stands out with effortlessness and elegance, with charm and caliber.



*The product images shown are for illustrative purpose only and not be an exact representation of the product. Accessories shown in the product image may not be part of the standard configuration.

We are interested in growing relationship



Technical Specifications

Overall Dimension*	2150 mm L x 1020 mm W
Platform Frame*	2080 mm L x 960 mm W
Mattress Platform*	1910 mm L x 865 mm W
Mattress Platform Height*	450mm to 830mm
Bed Extension*	215 mm
Backrest**	0° - 70°
Knee Rest**	0° - 45°
Trendelenburg**	0° - 17°
Reverse Trendelenburg**	0° - 17°
Safe Working Load	250 kgs
Patient Load Bearing Capacity^	200 kgs

Electrical Parameters

Power In (Mains Power Supply)	100 - 240 V ~, 50/60 Hz., Max 5A
Electrical Shock Protection	Class I
Degree of Shock Protection	Type B
Liquid Ingress Protection	IPX6
Duty Cycle	10 % (Two minutes for every eighteen minutes)

Safety Standards

General safety for medical beds	IEC 60601-2-52
General safety for medical electrical equipment	IEC 60601-1-1
EMI/EMC	IEC 60601-1-2

Variants


Features	EX6000 Plus with 380 & 392	EX6000 Plus with 360
Coverage span	Full length	¾ th length [§]

Both the variants are also available without Built-in Bed Extension feature.

* ±10 mm Engineering variation. ** ± 5° Engineering variation. ^When bed is configured with standard accessories.

[†]Factory fitted. [§]Can be converted into full length by adding 3rd set railing.

Add-On Accessories

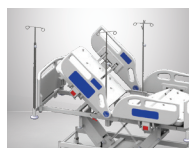

Heavy Duty and Telescopic IV Pole



O2 Cylinder Cage



Orthopaedic Traction Pulley



Stainless Steel Lifting Pole



Moulded Chart Holder



3rd Set of Safety Side Railings



Bed Extension[†]

- Single Section Premium Mattress and Four Section PU Foam Mattress
- Battery backup
- All variants are available with X-Ray permeable backrest (XR) feature as an add-on accessory, which are factory fitted.

Midmark (India) Private Limited

Head Office: Art Guild House, FOF, B-16, "B" Wing, 4th floor, Phoenix Market City, LBS Road, Kurla (West), Mumbai - 400 070
Toll Free: 1800 22 8020 | Tel: +91 22 4915 3000 | Fax: +91 22 4915 3100 | www.midmark.in | india@midmark.com

EX6000 Plus:12:07.01.24

Standard Features

- One touch key for cardiac chair, auto CPR, emergency head down
- Motorised backrest, kneerest, height adjustment, auto contour, trendelenburg & reverse trendelenburg
- Separate wired handset control for patient and attendant with lock out switches on attendant control panel
- Dual side manually operated CPR lever
- Calf elevation/DVT platform
- Angle indicator for backrest, upper leg, trendelenburg / reverse trendelenburg
- Polymer moulded split type safety side rails (2 sets)
- Polymer moulded head & foot panels, easily removable without using any tools
- Four corner safety rubber buffers
- 125 mm diameter twin wheel high grade synthetic body castors with central locking mechanism
- Inbuilt SMPS
- Easily removable polypropylene 5 section mattress top
- Built-in bed extension
- Urine bag holder
- Pull out linen holder
- Under bed night lamp

Provision for

- SS heavy duty IV pole at 4 corners
- Utility & telescopic IV pole above angle indicator

Accessories and color shown in the picture may not be part of standard configuration. Technical alterations reserved in view of our continuous improvements. Midmark India Pvt. Ltd. is the licensed user of the trademark "Midmark India" and the logo. Reproduction of this catalogue in any form is liable for legal action.