

OT Solutions



JANAK[®]
Because we care.



SURGISTAR⁺
MULTI PURPOSE ELECTROHYDRAULIC OPERATING TABLE



Provision for Ortho Attachment



User Friendly Handset



Railings with Locking System for Safety



Auxiliary Control on Column



Manual Override System Control



4 Heavy Duty Nylon Castors with Self-levelling Vacuum Assisted Manual Locking & Pump for Manual Override

Surgical procedures has become customizable with optimized workflow, improved efficacies yet cost-effective measures resulting in maximized profitability. Surgistar+ Electrohydraulic Operating Table systems are unmatched in versatility, with detachable and interchangeable head and leg sections for various surgical procedure, including radio translucent table top that offers outstanding C-arm imaging/fluoroscopy. Surgistar+ tables can be adapted to the desired operating environment and are capable of communicating with C-Arm systems to provide the ultimate in imaging and surgical performance. With unmatched adaptive positioning options, Surgistar+ facilitates seamless patient positioning and improves ergonomics for the surgical team.

Surgistar+ provides state-of-the-art solutions for cutting-edge, configurable surgical table systems. Surgistar+ multifunctional OT Tables enable flexible OR function for both minimally invasive and open surgical procedures.



Height Up & Down



Trendelenburg / Reverse Trendelenburg



Lateral Tilt (Left / Right)



Longitudinal Slide Towards Leg Side



Back section Adjustment

OT Solutions

Technical Specifications

Table Top Size (LxW)*	1930 x 533 mm
Minimum Height (Without Mattress)*	700 mm
Maximum Height (Without Mattress)*	1000mm
Trendelenburg / Reverse Trendelenburg**	25° / 25°
Lateral Tilt**	20° / 20°
Back Section Adjustment**	+75° / -25°
Longitudinal Table Top Slide (Optional)*	250mm
Head Section Adjustment**	+50° / -90°
Leg Section Adjustment**	+50° / -90°
Kidney Bridge*	125 mm
Weight Bearing Capacity^	200 kgs
Power Supply	220V AC - 50Hz

Standard Accessories

"L" Shaped Anaesthetic Screen	One pc
Padded Shoulder Support	One pair
Padded Lateral Support	One pair
Radio-Translucent Armboard with Pad	Two pcs
Goepel Knee Crutch with Straps	One pair
PU Foam Mattress	One set
Body Strap	One pc

Variant

Surgistar+ Core

4 Electrohydraulic Functions (Height, Trend/Rev Trend, Lateral Tilt and Back Section), Kidney Bridge, Single Leg Section, Handset, Auxillary Control on Column, Battery Backup, Manual Override System

Surgistar+ Pro

4 Electrohydraulic Functions (Height, Trend/Rev Trend, Lateral Tilt and Back Section), Kidney Bridge, Single Leg Section, Handset, Auxillary Control on Column, '0' Position Function, Flex-Reflex Position, Battery Backup, Manual Override System

Surgistar+ Advance

5 Electrohydraulic Functions (Height, Trend/Rev Trend, Lateral Tilt, Back Section and Longitudinal Slide of 250mm), Kidney Bridge, Single Leg Section, Handset, Auxillary Control on Column, Battery Backup, Manual Override System

Surgistar+ Max

5 Electrohydraulic Functions (Height, Trend/Rev Trend, Lateral Tilt, Back Section and Longitudinal Slide of 250mm), Kidney Bridge, Single Leg Section, Handset, Auxillary Control on Column, '0' Position Function, Flex-Reflex Position, Battery Backup, Manual Override System

Additional Accessories for Operating Table Surgistar+



For other accessories contact Midmark India sales representative

C-Arm systems to provide the ultimate in imaging and surgical performance




* ±10 mm Engineering variation. ** ±2° Engineering variation. ^When operating table is configured with standard accessories

Midmark (India) Private Limited

Head Office: Art Guild House, FOF, B-16, "B" Wing, 4th floor, Phoenix Market City, LBS Road, Kurla (West), Mumbai - 400 070
Toll Free: 1800 22 8020 | Tel: +91 22 4915 3000 | Fax: +91 22 4915 3100 | www.midmark.in | india@midmark.com

We are interested in growing relationships

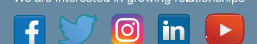

Surgistar+: 00:20:07.21