




AstirO

AstirO with SS fully collapsible safety side railing (pair)



01 X-Ray Permeable platform

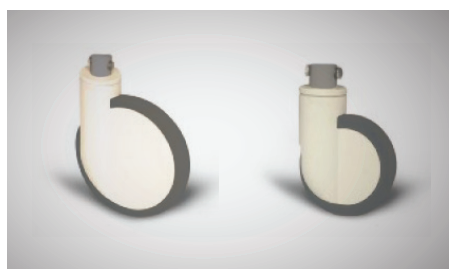

02 Multistage Switching Lock



01 Oxygen Cylinder Cage



01 Safety Head Board



01 Safety Side Railings



01 Safety corner buffers & push pull handles with rubber grippers (Available in red & blue colours)

Technical Specifications

Overall Dimensions*	2000 mm x 870 mm (L x W)
Stretcher Dimensions*	1840 mm x 660 mm (L x W)
Mattress Platform Dimensions*	1800 mm x 660 mm (L x W)
Backrest**	0° to 65°
Trendelenburg**	0° to 14°
Reverse Trendelenburg**	0° to 7°
Castor Diameter	200 mm / 125 mm
Corner Buffer Diameter	100 mm
Height* with 8" castor	710 mm to 910 mm
Height* with 5" castor	690 mm to 890 mm
Safe Working Load	200 Kgs
Patient Load Bearing Capacity^	175 Kgs

Standard Accessories

SS safety head board with parking bay
Oxygen cylinder cage
SS heavy duty telescopic IV pole with 2 hooks
Safety belts (pair)

Add on Accessories

Rexine covered two section PU foam mattress, 50 mm thick
Heavy duty IV pole (4 hooks)
Vertical oxygen cylinder cage (epoxy coated) with SS mounting rod (for small size cylinder weighing up to 10 Kgs)

AstirO is Engineered for Safe & Quick Intra Hospital Patient Transfer & Equipped for Emergency & Recovery Care Management

- Heavy duty non- staining, precision ball bearing swivel castors of 8" or 5" size with high-grade synthetic body for easy maneuverability
- Foot operated and centrally actuated lever at head and foot ends with multistage switching for total lock, direction lock or free mobility
- Gas spring assisted seamless adjustment of backrest, trendelenburg and reverse trendelenburg positions
- Convenient height adjustment on high quality hydraulic pump with foot operated levers on either sides of the trolley
- Two section mattress platform made of X-Ray permeable pre-laminated board, supported on strong tubular frame with removable stretcher
- Sliding multisize X-Ray cassette holder
- Individually operated SS safety side railings
- Safety corner buffers & push pull handles with rubber grippers (Available in red & blue colours)
- Provision of mounting heavy duty IV Pole at four locations
- Ample storage space

Safety Standards

Compliance for safety:
International Standards IEC 60601-2-52

Accessories and color shown in the picture may not be part of standard configuration. Technical alterations reserved in view of our continuous improvements. Midmark India Pvt. Ltd. is the Licensed user of the trademark "Midmark India" and the Logo. Reproduction of this catalogue in any form is Liabile for Legal action.

Midmark (India) Private Limited

Head Office: Art Guild House, FOF, B-16, "B" Wing, 4th floor, Phoenix Market City, LBS Road, Kurla (West), Mumbai - 400 070
Toll Free: 1800 22 8020 | Tel: +91 22 4915 3000 | Fax: +91 22 4915 3100 | www.midmark.in | india@midmark.com