

# Wardcare Solutions



**JANAK®**  
Because we care.

Fowler (Back and Knee / Leg adjustable) bed with conventional head and foot boards and side railing



Fowler (Back and Knee / Leg adjustable) bed with stainless steel laminated head and foot boards and side railing



Fowler (Back and Knee / Leg adjustable) bed with polymer moulded head and foot boards and side railing



Fowler (Back and Knee / Leg adjustable) bed with polymer moulded head and foot boards and side railing



## BX4000 Fowler Bed

The BX4000 fowler bed has the facility to raise the back and upper leg sections which are operated by the crank mechanism provided at the foot end of the bed. The raising section allows the nursing staff to change positions for various activities and clinical positions with relative ease.

\*The product images shown are for illustrative purpose only and not be an exact representation of the product. Accessories shown in the product image may not be part of the standard configuration.



# Wardcare Solutions

## Technical Specifications

Particulars	BX4000
Bed Frame Dimension*	2090 mm x 920 mm, 82" x 36" (L x W)
Mattress Platform*	1970 mm x 820 mm, 77" X 32" (L x W)
Buffer to Handle Dimension*	2310 mm 91"
Buffer to Buffer* (when castors provided)	2140 mm x 970 mm, 84" x 38" (L x W)
Backrest Angle**	0° - 70°
Knee Rest Angle**	0° - 50°
Safe Working Load	150 kgs
Patient Weight Bearing Capacity^	110 kgs

## Standard Features

- Raised backrest
- Raised kneerest
- Four section metallic perforated top
- Pre-treated and epoxy powder coated MS parts
- Provision for IV pole
- Removable handle
- Easily removable head and foot boards
- Supplied in SKD condition

## Safety Standards

General Safety for Medical Beds IEC 60601-2-52

## Add-On Accessories

### Railings and Legs

- Collapsible type epoxy powder coated railing (pair)
- Collapsible type epoxy powder coated MS safety side railing with SS tubes (pair)
- Swing down type epoxy powder coated MS safety side railing (pair) (factory fitted)
- Collapsible type epoxy powder coated vertical safety side railing (pair)
- Collapsible type full-length type SS railing with epoxy powder coated frame (pair)
- Single set big size polymer railing (factory fitted)
- Split type polymer moulded safety side railing (pair) (factory fitted)
- H' type legs 910 mm (36") width & 600 mm (24") height with rubber shoes
- H' type legs 910 mm (36") width & 450 mm (18") height with rubber shoes
- H' type legs 910 mm (36") width & 600 mm (24") height with 125 mm high grade synthetic castors, two with brake and two without brake
- H' type legs 910 mm (36") width & 450 mm (18") height with 125 mm high grade synthetic castors, two with brake and two without brake

- Four section mattress with 4" thick PU foam of 40 density covered with PVC
- Telescopic IV pole with 2 hooks
- Telescopic IV pole with 4 hooks
- Orthopaedic traction pulley for trauma cases
- Fixed polymer moulded handle (set of 2)
- Removable handle
- Moulded chart holder
- Mattress arrestor (set of 2)
- Mosquito net poles
- 4" safety rubber buffer at four corners
- Urine bag holder

### Head and Foot Boards

- Polymer moulded head and foot boards
- SS laminated head and foot boards
- Conventional head and foot boards


\*± 10 mm Engineering Variation and \*\* ± 5° Engineering Variation ^ When Bed is configured with standard accessories  
Note: Option of fixed polymer moulded handle available in lieu of removable handle.



Accessories and color shown in the picture may not be part of standard configuration. Technical alterations reserved in view of our continuous improvements. Midmark India Pvt. Ltd. is the licensed user of the trademark "Midmark India" and the logo. Reproduction of this catalogue in any form is liable for legal action.

## Midmark (India) Private Limited

Head Office: Art Guild House, FOF, B-16, "B" Wing, 4th floor, Phoenix Market City, LBS Road, Kurla (West), Mumbai - 400 070  
Toll Free: 1800 22 8020 | Tel: +91 22 4915 3000 | Fax: +91 22 4915 3100 | www.midmark.in | india@midmark.com

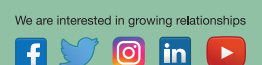

We are interested in growing relationships

BX4000:09:20:09.22